


DROPLET + CONTACT ISOLATION

DONNING- SEQUENCE FOR PUTTING ON PERSONAL PROTECTIVE EQUIPMENT (PPE)

Perform Donning **in the Hallway**, prior to entering the room.


1. HAND HYGIENE:

- Use hand sanitizer to clean hands

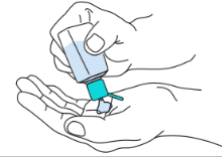

2. GOWN:

- Fully cover torso from neck to knees, arms to end of wrists, and wrap around the back
- Fasten ties at back of gown

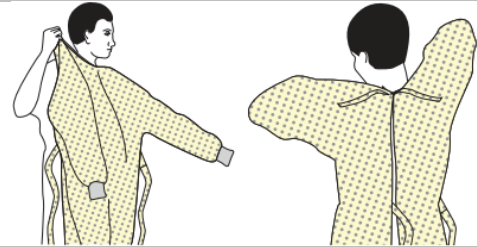

3. SURGICAL MASK:

- Secure tie or elastic bands at middle of head and neck
- Fit flexible band to nose bridge
- Fit snug to face and below chin
- Mask with attached face shield is acceptable to use (not pictured); if using, no additional eye protection required




4. EYE PROTECTION (GOGGLES OR FACE SHIELD):

- Place over face and eyes and adjust to fit




5. HAND HYGIENE:

- Use hand sanitizer to clean hands

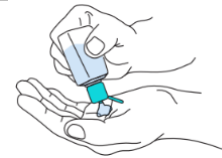

6. GLOVES:

- Extend to cover wrist of isolation gown

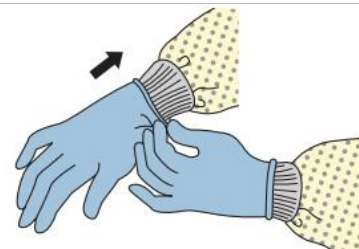

USE SAFE WORK PRACTICES TO PROTECT YOURSELF AND LIMIT THE SPREAD OF CONTAMINATION

- Keep hands away from face
- Limit surfaces touched
- Change gloves when torn or heavily contaminated
- Perform hand hygiene

ADAPTED FROM:




DROPLET + CONTACT ISOLATION

DOFFING- HOW TO SAFELY REMOVE PERSONAL PROTECTIVE EQUIPMENT (PPE)

- Remove Gown and Gloves in the *Patient Room*.
- Remove Eye Protection and Mask in the *Hallway*.
- Safe removal of PPE is essential to preventing contamination of clothing, skin, or mucous membranes with potentially infectious materials

Remove PPE in the following sequence:**1. In Room: REMOVE GOWN AND GLOVES:**

- Gown front and sleeves and the outside of gloves are contaminated!
- If your hands get contaminated during gown or glove removal, immediately wash your hands or use an alcohol-based hand sanitizer
- Grasp gown in front and pull away from your body so that the ties break, touching outside of gown only with gloved hands
- While removing the gown, fold or roll the gown inside-out into a bundle
- As you are removing the gown, peel off your gloves at the same time, only touching the inside of the gloves and gown with your bare hands
- Place the gown and gloves into a waste container

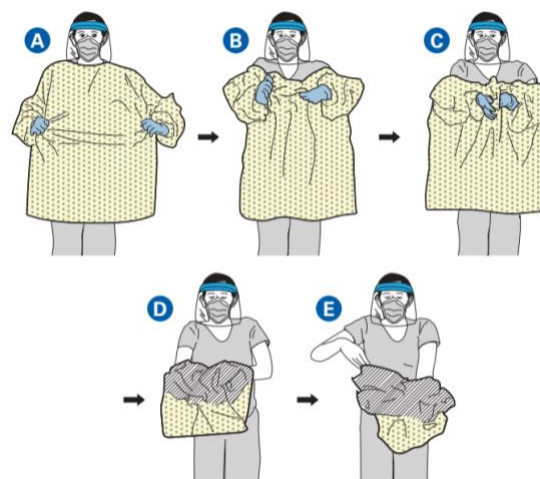**2. In Room: HAND HYGIENE:**

- Use hand sanitizer to clean hands prior to removing remaining PPE

**3. In Hall: Remove EYE PROTECTION (GOGGLES OR FACE SHIELD):**


- Outside of goggles or face shield are contaminated!
- If using mask with attached face shield, skip to "5. Surgical Mask" doffing below
- Remove goggles or face shield from the back by lifting head band and without touching the front of the goggles or face shield
- Discard in a waste container

**4. In Hall: HAND HYGIENE:**

- Use hand sanitizer to clean hands prior to removing remaining PPE

**5. In Hall: Remove SURGICAL MASK:**

- Front of mask is contaminated — DO NOT TOUCH!
- If your hands get contaminated during mask removal, immediately perform hand hygiene
- Release bottom ties or elastics of the mask, then the ones at the top, and remove without touching the front
- Discard in a waste container

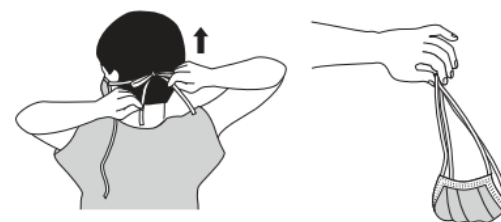**6. In Hall: WASH HANDS OR USE ALCOHOL-BASED HAND SANITIZER**